

Abutment Kit **AK-1**

Mill Custom Titanium Abutments

Save Time. Save Money. Save Outsourcing.

Available for the DWX-42W




DWX-42W



With a better fit and greater overall esthetics, custom titanium abutments are fast becoming the standard choice for clinicians and dental laboratories. The AK-1 Abutment kit for the DGSHAPE DWX-42W wet dental mill brings the production and milling of custom titanium abutments in-house. It eliminates the need to outsource and provides you with a quality, full-service implant restoration solution.


Combine a DWX-42W with the AK-1 Abutment Kit that features selected equipment for use with your choice of the industry's most popular attachments and pre-milled titanium blanks*. This out-of-box solution provides all the milling burs, filters and tools you need to set-up your DWX-42W for custom titanium abutments.

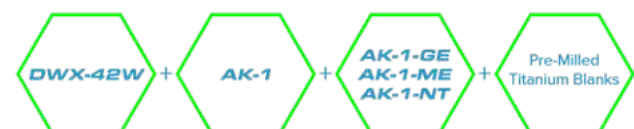

*AK-1 material adapter attachments and pre-milled titanium blanks sold separately.

Custom Titanium Abutments

Using a kit that includes specially developed accessories and tools, mill custom abutments in just 3-easy-steps with the DWX-42W.

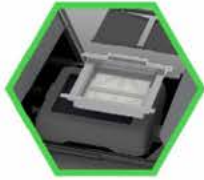

1

Switch-Out Tray

Replace the regular tank with the special AK-1 titanium tray.



2

Connect AK-1 Material Adapter

Connect your choice of AK-1 material adapter* to the work material holder.



3

Mill & Complete

Insert and attach the pre-milled titanium block, output your milling job and complete.

*Material adapter sold separately

Your AK-1 Abutment Kit Includes



Three Types of Special-Purpose Milling Burs

Specially engineered for milling GeoMedi, MEDENTIKA, and NT Trading pre-milled blocks.



Switch-Out Tray and Special Purpose Filters

Designed specifically for the collection of titanium milling waste—supplied with two replacement filters.



Torque Driver and Hexagonal Wrench

Fix materials to the adapter with the appropriate levels of torque and make other adjustments.



Receptacle Tray

Keep your milling materials close-to-hand and organized.

Material Adapter Attachments

(Sold Separately)



AK-1-GE

Compatible with GeoMedi Geo CMFit® series



AK-1-ME

Compatible with MEDENTIKA PreFace® series



AK-1-NT

Compatible with NT-Trading NT-PREFORM® series

Highlights

- Low investment and low cost per unit combines with in-house production savings for an exceptionally economical solution.
- Clinicians and lab technicians have full control over design and materials from start-to-finish—mill both the custom abutment and implant supported restoration in-house.
- Custom Titanium Abutments are more versatile—offering flexibility for a passive fit and prosthetic correction of implant angulation, as well as greater esthetics, especially in the anterior region.
- Save time while providing a better experience for dental patients—custom abutments require fewer chairside visits and adjustments.
- A simple update to the latest version of MillBox Software for your DWX-42W offers you an intuitive CAM solution for custom abutment milling.
- Production with your chosen AK-1 material adapter and pre-milled titanium blanks on the DWX-42W is compatible with some of today's most popular CAD software systems.



Authorized Premier Distributor

ETI DIGITAL TECHNOLOGY

Complete Source for the Future

Direct (714) 238-1490

www.etidigitaltechnology.com

sales@etidigitaltechnology.com



FOLLOW US!

COMPLETE SUPPORT FROM OUR EXPERIENCED DIGITAL TEAM